

EtCO₂ Accessories

INFINIUM[®]



Emergency & Trauma

Sedation & Anesthesia

Critical Care

Lung Function

Modules . Cannula . Sampling Lines . Filters

Philips Respironics EtCO2 Modules

Philips Respironics Capnostat 5 Mainstream EtCO2

Infinium Part: 152.001.0002

Respironics Part: M2501A

Uses sophisticated infrared absorption spectroscopy to measure EtCO₂. Measurements are taken at the patient's airway, so response is faster and there is less chance of erroneous artifact data. With no moving parts the sensor is durable and provides reliable readings practitioners can count on.

- Direct airway measurement provides response times less than 60ms
- Easy to use with simple snap-in sensors. No moving parts
- Very durable and can withstand multiple drops
- No regular calibration required



Philips Respironics LoFlo Sidestream EtCO2 Module

Infinium Part: 152.001.0001

Respironics Part: M2741A

Uses sophisticated infrared absorption spectroscopy to measure EtCO₂ & respiration rate. The sensor delivers a precise capnogram for respiratory rates up to 150 breaths per minute. With a sampling rate of 50 ml/min, it is ideal for use with infants and smaller patients.

- Waveform data in less than 20 seconds
- Easy to use with simple snap-in sensors. No moving parts
- Very rugged, shock resistant and IPX4-splash-proof
- No regular calibration required





Respironics EtCO2 Parts & Accessories

	Description	Infinium Part	Qty	Respironics Reference
	Respironics Adult Straight Sample Line with Adapter	300.001.0002	BX/10	3472ADU, 989803206711, 989803144521, M2768A
	Respironics Straight Sample Line with Humidifier	300.001.0003	BX/10	3475, 98980314456, M2776A1
	Respironics Adult Nasal Sample Line	300.001.0004	BX/10	3468ADU, 989803206611, 989803144421, M2744A
	Respironics Mainstream Airway Adapter Adult	300.001.0005	BX/10	6063, 989803142661, M2533A
	Respironics Mainstream Airway Adapter Neonate	300.001.0006	BX/10	6312, 989805618811
	Respironics Adult Straight Sample Line with Humidifier and Adapter	300.001.0007	BX/10	3473ADU, 989803206721, 989803144531

Luer Lock CO2 & Oxygen Cannula Sets



		Description	Infinium Part	Qty
		CO2 Nasal Cannula with (male connector)	300.001.0010	BX/25
		CO2 Nasal Cannula with (female connector)	300.001.0011	BX/25
		Co2 Nasal Cannula with O2 tube (male connector)	300.001.0012	BX/25
		Co2 Nasal Cannula with O2 tube (female connector)	300.001.0013	BX/25
		1-inch CO2 sampling line filter	300.001.0016	BX/25



Oridion - Capnoline EtCO2 Sampling Lines

Comparable Brands / Styles: Oridion, Microstream, Capnoline

Compatible With: Medtronic Capnostream with Microstream - Philips Intellivue & Efficia Series with Microstream Capnography

Adult Intubated EtCO2 Sampling Line for Short Term Use

Infinium Part: 300.010.002
Microstream Ref: MVAI
Tubing: STD 6.5ft
Packaging: EACH

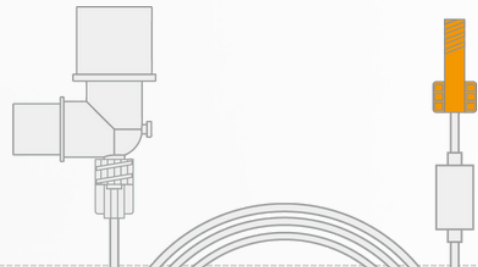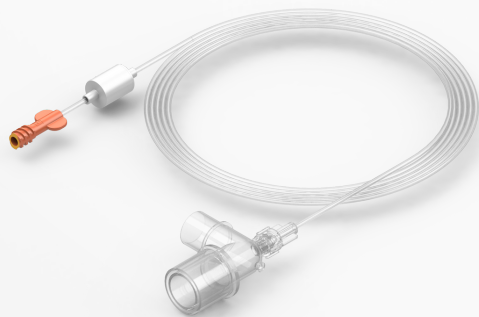

Adult Nasal Sampling Line for Short Term Use

Infinium Part: 300.010.001
Microstream Ref: MVAN
Tubing: STD 6.5ft
Packaging: EACH



Adult Nasal Sampling with O2 Delivery for Short Term Use

Infinium Part: 300.010.003
Microstream Ref: MVANO
Tubing: STD 6.5ft
Packaging: EACH






Infinium EtCO₂ Monitoring Modules & Accessories

With over 30 years of patient monitoring manufacturing, Infinium understands the importance of high-quality and cost effective capnography solutions. Talk with us about how to improve your daily operations.



INFINIUM®

Infinium Medical
12151 62nd St North #5
Largo, FL 33773 • USA

 (1) 866-918-8434
 infiniummedical.com
 sale@infiniummedical.com